

JORDAN AND HAMBURG LLP
122 East 42nd Street
New York, NY 10168
Tel: (212) 986-2340
Fax: (212) 953-7733

JC03 Rec'd PCT/PTO 26 APR 2001

09/830511

Docket No. F-6954
Date April 26, 2001

"EXPRESS MAIL" MAILING LABEL NO. EL 666 237 967 US

DATE OF DEPOSIT April 26, 2001

I hereby certify that this paper or fee is being deposited with the United States Postal Service "Express Mail Post Office to Addressee" service under 37 CFR 1.10 on the date indicated above and is addressed to the Assistant Commissioner for Patents, Washington, D. C. 20231.

Madeline Gonzalez

(Typed or printed name of person mailing paper or fee)

Madeline Gonzalez
(Signature of person mailing paper or fee)

THE ASSISTANT COMMISSIONER FOR PATENTS
Washington, D. C. 20231

[] ATTN: BOX PATENT APPLICATION

[X] ATTN: BOX PCT

[X] THIS IS THE NATIONAL STAGE OF PCT/EP99/08178 FILED October 28, 1999

Sir:

Transmitted herewith for filing is the [X] Utility [] Design patent application of:

Inventor/Application Identifier: Alexandre TRETIKOV et al.

For: **ULTRATHIN-WALLED MULTIWELL PLATE FOR
HEAT BLOCK THERMOCYCLING**

Enclosed are:

- [X] 2 sheets of drawings ([X] formal [] informal size A4).
[X] 9 pages of specification, including claims and abstract.
[X] 11 total pages
[X] Combined Declaration/Power of Attorney
[X] Newly executed
[] Copy from prior application
[] Inventors deleted; see attached statement
[] Inventor Information Sheet
[] Incorporation By Reference. The entire disclosure of the prior application, from which a copy of the oath or declaration is supplied, is considered as being part of the disclosure of the accompanying application and is hereby incorporated by reference therein
[] Sequence Listing
[] Computer Readable Copy
[] Paper copy
[] The undersigned hereby affirms that the content of the paper and computer readable copies of the Sequence Listing are the same.
[] Cancel in this application original claims _____ of the prior application before calculating the filing fee.

CLAIMS FILED

For	Number Filed	Number Extra	Rate	Basic Fee \$860.00
Total Claims	7	(over 20)	x \$18.00	
Independent Claims	1	(over 3)	x \$80.00	
[] Multiple Dependent Claim			\$270.00	
[] Reduce by 50% for Small Entity			-	
[] Foreign Language Filing Fee			\$130.00	

TOTAL FILING FEE **\$860.00**

[X] Please charge Deposit Account No. 10-1250 in the amount of
A duplicate copy of this sheet is attached.

\$860.00.

[X] Please charge to Deposit Account No. 10-1250 any further fees under
37 CFR 1.16; 37 CFR 1.17; 37 CFR 1.492.

☒ Return Receipt Postcard

☒ Preliminary Amendment

☒ Assignment to Hans-Knoell-Institut fuer Naturstoff-Forschung e.V.

☐ Assignment is of record in prior application Serial No. _____

☒ Assignment Recordation Form Cover Sheet.

☒ Charge \$40.00 to Deposit Account No. 10-1250 for recording Assignment.

☒ Information Disclosure Statement and/or Information Disclosure Citation

☐ English translation

☐ Small Entity Status is asserted.

☐ Applicant hereby claims the benefit of the filing date of the following provisional application(s) under the provisions of 35 USC 119.


☒ Applicant hereby claims the benefit of the filing date of the following applications under the provisions of 35 USC 119 of which certified copies ☐ will follow ☐ are enclosed ☒ have been filed in the International Bureau ☐ were filed in prior application No. _____

Europe Patent Appln. No. 98120187.4 filed October 29, 1998.

☐ This is a ☐ Continuation ☐ Divisional ☐ Continuation-in-Part of prior application Serial No. _____

☐ Amend the specification by inserting before the first line the sentence:
--This is a ☐ continuation, ☐ division, ☐ continuation-in-part, of application
Serial No. _____, filed _____ --

JORDAN AND HAMBURG LLP

By 
C. Bruce Hamburg
Reg. No. 22,389
Attorney for Applicants

09/830511

JC16 Rec'd PCT/PTO 26 APR 2001

F-6954

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicant : Alexandre TRETIKOV et al.
Serial No. : Not yet known (U.S. National Stage of
PCT/EP99/08178 filed October 28, 1999)
For : ULTRATHIN-WALLED MULTIWELL PLATE
FOR HEAT BLOCK THERMOCYCLING

Assistant Commissioner for Patents
Washington, D.C. 20231

PRELIMINARY AMENDMENT

Sir:

Preliminary to examination, please amend this application as follows:

IN THE CLAIMS:

Cancel claim 6.

Amend claim 1 as follows:

1. (Amended) Ultrathin-walled multiwell plate for heat block thermocycling of samples comprising an array of small-volume wells of identical height with the similarly shaped sample wells formed in the top surface of the heat block of the thermocycler, wherein the walls of the wells have an average thickness of 20-40 microns.

The Appendix hereto shows the changes in claim 1.

09830511-042301
T03240 T5000050


F-6954

REMARKS

This corresponds to an amendment filed in the international stage.

Respectfully submitted,

Jordan and Hamburg LLP

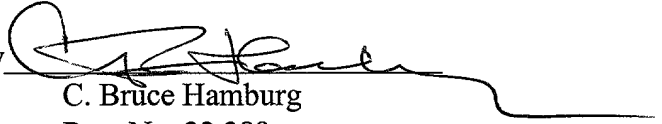By 
C. Bruce Hamburg
Reg. No. 22,389
Attorney for Applicants

122 East 42nd Street
New York, New York 10168
(212) 986-2340

1093054-044304
1093054-044304

APPENDIX

1. (Amended) Ultrathin-walled multiwell plate for heat block thermocycling of samples comprising an array of small-volume wells of identical height with the similarly shaped sample wells formed in the top surface of the heat block of the thermocycler, wherein the walls of the wells have an average thickness of 20-40 microns.

09030301.043601

Ultrathin-walled multiwell plate for heat block thermocycling

The invention relates to plastic plates for conventional heat block thermocycling
5 of biological samples, particularly to multiwell plates. More specifically, it relates to
ultrathin-walled multiwell plates with an improved heat transfer to small-volume
samples. Such plates can be used for rapid temperature cycling of multiple, small-volume
samples (i.e. 1-20 μ l) by using heat block thermocyclers with an increased block
temperature ramping rate (i.e. 4° C/second and greater) and standard heated-lid
10 technology for sealing the samples.

Temperature cycling of biological samples is a central moment in DNA
amplification by the polymerase chain reaction (PCR) (Saiki et al., Science, 239, 487-491
[1988]). Much effort is being expended in developing various alternative reactors and
technologies for rapid temperature cycling of small-volume samples (Kopp et al., Science
15 280, 1046-1048 [1998]; Belgrader et al., J. Forensic Science 43, 315-319 [1998]; Wittwer
et al., Analytical Biochem., 186, 328-331 [1990] and U.S. Patent No 5,455,175; Woolley
et al., Analytical Chem., 68, 4081-4086 [(1996)]).

One commercially available type of microreactor and thermocycler for rapid
temperature cycling of small samples is a glass capillary tube and a hot-air thermocycler
20 from Roche Molecular Biochemicals (cat No. 1909 339 and cat No. 2011468, respectively).
The glass capillary tube can hold reaction volumes ranging from 10 to 20 μ l. The hot-air
thermocycler can hold 32 capillaries and perform 30 - 40 PCR cycles in 20-30 minutes.
However, these rapid DNA amplification technology is connected with various
disadvantages, for example:

- 25 a) The handling of the individual capillaries is relatively cumbersome.
- b) The relatively large glass surface adsorbs components of the standard PCR-mixtures.
This might inactivate the reaction. Therefore, various carrier molecules, i.e. proteins
or even DNA, must be added and the concentrations of the components reoptimized.
- c) The cost of the capillary tube, as a disposable PCR container, is high when compared
30 to the standard 0.2 ml PCR tube.
- d) The experimental throughput using this system is limited.

It is surprising that only little research has been conducted to improve the basic
performance in sample size and speed of the widely used, conventional heat block

thermocycling of samples contained in plastic tubes or multiwell plates. One known improvement of heat block temperature cycling of samples contained in plastic tubes has been described by Half et al. (Biotechniques, 10, 106-112, [1991] and U.S. Patent No 5,475,610). They describe a special PCR reaction-compatible one-piece plastic microcentrifuge tube, i.e. a thin-walled PCR tube. The tube has a cylindrically shaped upper wall section, a relatively thin (i.e. approximately 0.3 mm) conically-shaped lower wall section and a dome-shaped bottom. The samples as small as 20 µl are placed into the tubes, the tubes are closed by deformable, gas-tight caps and positioned into similarly shaped conical wells machined in the body of the heat block. The heated cover compresses each cap and forces each tube down firmly into its own well. The heated platen (i.e. heated lid) serves several goals by supplying the appropriate pressure to the caps of the tubes: it maintains the conically shaped walls in close thermal contact with the body of the block; it prevents the opening of the caps by increased air pressure arising in the tubes at elevated temperatures. In addition, it maintains the parts of the tubes that project above the top surface of the block at 95° -100° C in order to prevent water condensation and sample loss in the course of thermocycling. This made it possible to exclude the placing of mineral oil or glycerol into the wells of the block in order to improve the heat transfer to the tubes and the overlaying of the samples by mineral oil that prevented evaporation but also served as added thermal mass. In addition, the PCR tubes can be put in a two-piece holder (US patent 5,710,381) of an 8x12, 96-well microplate format, which can be used to support the high sample throughput needs with any number between 1 and 96 individual reaction tubes.

In DE 4022792 the inventors describe a plate with cylindrically shaped walls of the wells and spherically shaped bottoms thereof. The individual wells of the plate were formed by melting a polycarbonate sheet in the range of 0.27-0.5 mm by a stream of hot air. This technology leads to relatively thin walls in the range of 0.08-0.2mm. The biological samples were placed into the wells, covered with polycarbonate film (0.1 mm) and the individual wells were thermosealed by a special press. Upon sealing the plate was placed on the thermoblock and fixed by screws. Though theoretically the heat transfer to the samples is improved, however, the way of positioning the plate on the block and the cylindrical and spherical geometry of the well prevent a close thermal contact with the heating block. During thermocycling, due to the large thermal expansion, the plate fixed by

screws becomes deformed and the close thermal contact is not maintained anymore.

Therefore, by using the above technology rapid cycling reactions cannot be performed.

The other known improvement of heat block thermocycling is described in PCT patent application WO 98/43740. It concerns a heat block thermocycler with an increased
5 ramping rate, i.e. $4^{\circ}\text{C}/\text{second}$). The thermocycler can hold 96 PCR tubes (each of a volume of 0.2 ml) or 96-well PCR plates. Theoretically, the thermocycler can perform 30 PCR cycles in 20-30 minutes, provided that only a few seconds are spent to reach the temperature equilibrium between the heat block and the samples.

However, as described in U.S. Patent No 5,508,197, even if the temperature of the
10 heat-transfer media, i.e. water, is changed almost instantaneously, it takes approximately 15 seconds to reach equilibrium between water and the 15-20 μl samples in the standard PCR plates. This means that for 30 PCR cycles approximately 20 minutes are spent to reach the equilibrium between heat-transfer media and the 15-20 μl samples in the plates.

In comparison, the above mentioned heat block cycler (WO 98/43740) operating
15 at a ramping rate of $4^{\circ}\text{C}/\text{second}$, needs for the heat-block temperature transitions during 30 PCR cycles 10 minutes only. This shows that the major limiting factor for rapid temperature cycling of small samples in plastic PCR tubes or PCR plates is the low efficiency of the heat transfer through the walls of conventional PCR tubes or plates, respectively.

20 The present invention concerns plastic multiwell plates for performing heat block thermocycling of multiple samples. More specifically, it concerns ultrathin-walled multiwell plates with an improved heat transfer to small samples. Ultrathin-walled multiwell plates are suited for rapid, oil-free, heat block temperature cycling of small-volume samples (i.e. approximately 1-20 μl), whereas the lower limit is given by the
25 reliability of the conventional pipetting systems.

Figure 1 illustrates an example of a multiwell plate according to the invention.

Figure 2 illustrates the positioning of the plate in the block of the thermal cycler.

One aspect of the present invention concerns the considerably decreased thickness (i.e. approximately 7.5-15 fold) of the well walls when compared to known thin-walled
30 PCR tubes (U.S. Patent No 5,475,610). This can be reached, for example, by means of thermoforming the plates out of thin thermoplastic films. Such thermoplastic films are, for example, polyolefin films, such as metallocene-catalyzed polyolefin films and/or copolymer films. Usually, the multiwell plate is vacuumformed out of cast, unoriented

09330541.042601
FO9240" T509660

polypropylene film, polypropylene-polyethylene copolymer films or metallocene-catalyzed polypropylene films. The film is formed into a negative ("female") mould comprising a plurality of spaced-apart, conically shaped wells which are machined in the body of a mould in the shape of rectangular- or square-array. The thickness of the film
 5 for vacuumforming conically shaped wells is chosen according to the standard rule used for thermoforming, i.e. thickness of the film = well draw ratio x thickness of the wall of the formed well.

For example, vacuumforming wells with a draw ratio of two and an average thickness of the walls of 30 microns results in a film thickness of 60 microns. The
 10 average optimum wall thickness was found to be 20-40 microns. The thickness of the well is reduced 7.5-15 fold when compared to the wall thickness of the formerly improved PCR tube described in U.S. Patent No 5,475,610. Using the Fourier equation for heat transfer and the equation for temperature transfer through solid substances, it can be shown that heat transfer through one square millimeter of the surface of the well of the
 15 plate is increased 7.5-15 fold and the time of temperature transfer through the wall is decreased 56-225 fold when compared to the said PCR tube. This drastic decrease in time can be explained by the fact that the time needed for the transfer of temperature front is proportional to the square power of distance. It can be easily calculated that the time of the temperature transfer through the ultrathin walls of the multi-well plate is in the range
 20 of milliseconds, whereas for the said PCR tube (U.S. Patent No 5,475,610) it is in the range of seconds. This explains the well known fact that thin (20-40 microns) plastic films are poor thermo insulators.

The thickness of the walls of the formed wells is gradually reduced to the bottoms of the wells due to vacuumforming of the wells into a negative mould. This geometry of
 25 the walls of the wells provides several advantages:

- The relatively thick upper parts of the walls of the wells cause additional rigidity of the whole multiwell plate.
- During heating of the heat block of the thermocycler, a vertical temperature gradient is formed in the sample, due to the gradient of the well-wall thickness. This vertical
 30 temperature gradient causes intensive convective mixing of the sample in conically shaped wells and increases the heat transfer through the sample. In comparison, this convective mixing of the sample is much less efficient in conventional PCR plates/tubes with a uniform wall thickness.

Another aspect of the invention concerns the height of the wells of the multiwell plate. The height of the conically shaped wells is equal to the height of the similarly shaped sample wells machined in the body of the heat block. Thus, this geometry of the wells (2) enables the positioning of the plate (1) on the heat block (4) as shown in Figure 2. As shown (Figure 2), in contrast to the conventional PCR plates, the walls of the wells (2) of the multi-well plate (1) do not project above the top surface of the block (4). The type of positioning provides several advantages: The pressure caused by the screw (12) to the lid (10) (heating element (11)) can be increased in order to obtain efficient sealing of the samples (9) sealed, for example, by a silicon mat (13). In this case the pressure is actually directed to those parts of the multiwell plate (1) which are supported by the top surface of the heat block (4) (or by parts of the top surface surrounding individual wells depending on the geometry of the heat block) and not to the thin walls of the wells of the plate as it is the case for the PCR tubes or conventional PCR plates. This advantage makes it possible to increase the sealing pressure of the heated lid (10) several fold when compared to the conventionally used pressure of 30-50 g per well without cracking the conically shaped walls of the wells (2).

The extremely thin walls of the wells, i.e. 20-40 microns, are highly flexible as the multiwell plates are thermoformed out of highly elastic films (or sheets depending on the draw ratio). The walls of the wells are highly resistant against stress cracking, due to their flexibility and elasticity. As the wells of the plate, positioned on the heat block, are tightly sealed at room temperature, the air pressure in the wells will increase at elevated temperatures. The increased air pressure causes a deformation of the walls of wells and brings them in tight thermal contact with the surface of the walls of the individual sample wells machined in the body of the heat block. Standard PCR plates (having relatively thick and rigid walls of the wells) require that the conically shaped walls of the wells have to match perfectly with the shape of the wells machined in the body of the heat block to guarantee a close thermal contact (see for example U.S. Patent No 5,475,610). This requirement is not as critical for the ultrathin walled multiwell plates of the invention, due to flexibility and elasticity of the walls of the wells. Using this advantage, special shapes of both, the walls of the wells of the plate and the wells of the heat block can be differently designed. These differently designed wells can promote an even closer thermal contact after positioning the plate into the heat block.

Another aspect of the invention concerns the frame of the multiwell plates. As the plates can be formed of very thin films (depending on the draw ratio of the well; supra) the flexibility of, for example, standard-format plates, i.e. 96-well PCR (8,5 x 12,5 cm) plates, is such that handling is not easily possible anymore. Therefore, depending on the geometry of the plate, a supporting frame might be needed, for example for industry standard formats, i.e. 96-, 192-, 384-well PCR plates. This frame can support, for example in case of small plates, the edges of the plate, or individual wells of the plate, or groups of wells. For handling with robots, for example, the frame can be injection molded in the form of the standard skirted microplates containing the array of holes in the top surface of the frame matching the array of wells of the ultrathin multiwell plate. The plate can be attached to the frame by for example heat bonding. However, for small format plates including the frame can be formed as a single piece by using specially designed moulds.

The polypropylene-based plastics are PCR-compatible and therefore widely used for injection molding of PCR tubes and/or multiwell plates. In addition, they are resistant to stress cracking and have a reduced water vapor sorption when compared to other plastics (e.g. polycarbonate). Such plates can be thermoformed in both, standard industry formats, i.e. 96-, 192- and 384-well PCR plates for large scale applications, supported by robots and small foot-print formats to match small foot-print thermocyclers, i.e. "personal thermocyclers".

The following example serves to illustrate the invention but should not be construed as a limitation thereof.

25 Example:

Fig.1 illustrates a 36-well ultrathin walled multiwell plate according to the invention. The plate was designed for rapid temperature cycling of samples ranging from 0.5-4 μ l using a small foot-print peltier-driven heat block thermocycler supplied with a "wine-press" type heated lid (Fig. 2). The volume of the wells is 16 μ l and the distance between the wells is 4.5 mm, i.e. industry standard for high sample density 384-well PCR plates. The diameter of the openings of the wells is 3.8 mm and the height of the wells is 3 mm. The average thickness of the walls of the wells is 30 μ m. The frame (3) was cut out of a polypropylene sheet of a thickness of 0.5 mm and heat bonded to the plate (1). The area of the plate (1)

is 30 x 30 mm. As shown in Figure 1, the handling of the plate (1) containing the multiple wells (2) is facilitated, by a rigid 0.5-1 mm thick plastic frame (3) which is heat bonded to the plate. As shown in Figure 2, the frame (3) is not in direct thermal contact with the block (4) during thermocycling because the inner contour (5) of the frame (3) matches the outer contour (6) of the heat block (4) of the thermocycler (7 = thermoelectric heat pump and 8 = air-forced heat sink).

The ultrathin walled multiwell plate according to the invention (Fig.1) was experimentally tested for the amplification of a 455-base pairs long fragment of human papilloma virus DNA. The sample volume was 3 μ l. For various PCR reactions, the average ramping rate of the thermo cycler was varied from 4° C to 8° C per second. The samples (i.e. standard PCR-mixtures without any carrier molecules) were transferred into the wells of the plate by means of conventional pipetting equipment. The plate was covered by standard sealing film (Microseal A; MJ-Research, USA), transferred into the heatblock of the thermocycler and tightly sealed by the heated lid as shown in Fig. 2. Upon sealing, a number of 30 PCR cycles was performed in 15-25 minutes depending on the ramping rate of the thermo cycler. The PCR product was analyzed by conventional agarose electrophoresis. The 455-base pairs long DNA fragment was amplified with a high specificity at the indicated ramping rates (supra).

Plates according to the invention with well volumes of 35 μ l were successfully tested for temperature cycling of samples of a volume of 20 μ l. Thereby, 30 PCR cycles were performed in 20-30 minutes at a ramping rate of 6° C per second. Surprisingly, although the average thickness of the walls was 20 microns and the volume of the wells was 35 μ l, samples of a volume of as few as 0.5 μ l can be easily amplified without reducing the PCR efficiency.

In conclusion, the ultrathin walled multiwell plates according to the invention, allow a simple and rapid loading of multiple samples by conventional pipettes, rapid sealing of all samples by using conventional sealing films and rapid DNA amplification (15-30 minutes for 30 cycles) with an improved specificity typical for rapid cycling (Wittwer et al., Analytical Biochem., 186, 328-331 [1990]) using appropriate heat block thermocyclers (i.e. ramping rate in the range of 4° C to 8° C per second).

MT 34 AMEND

Claims

1. Ultrathin-walled multiwell plate for heat block thermocycling of samples comprising an array of small-volume wells of identical height with the similarly shaped sample wells formed in the top surface of the heat block of the thermocycler, wherein the walls of the wells have an average thickness of 20-40 microns.
2. Ultrathin-walled multiwell plate according to claim 1, wherein the height of the wells of the plate is not more than the height of the sample wells formed in the top surface of the heat block of the thermocycler
3. Ultrathin-walled multiwell plate according to claim 1, wherein the walls of the wells are conically shaped.
4. Ultrathin-walled multiwell plate according to claim 1, wherein the thickness of the walls of the wells decreases from top to bottom.
5. Ultrathin-walled multiwell plate according to claim 1, wherein the wells of said multiwell plate are thermoformed into negative mould.
6. Ultrathin-walled multiwell plate according to claim 1, wherein the walls of the wells are deformable.
7. Ultrathin-walled multiwell plate according to claim 1, wherein the said microwell plate comprises a rigid supporting frame.
8. Ultrathin-walled multiwell plate according to claim 1, wherein the volume of the well is in the range of 16-85 μ l.

09930511.042601
T09240" T50E360

Ultrathin-walled multiwell plate for heat block thermocycling

Abstract of the disclosure

Ultrathin-walled multiwell reactors for heat block thermocycling of samples comprising an array of small-volume wells of identical height with similarly shaped sample wells formed in the top surface of the heat block of the thermocycler are provided. The multiwell plates are preferentially vacuumformed out of a 30-50 micron thick thermoplastic film and can be used for rapid, oil-free temperature cycling of small (1-10 μ l) volume samples.

09830511 042601

Figure 1



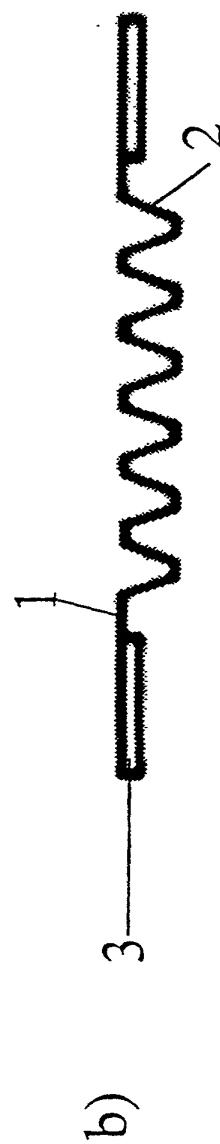


Figure 2

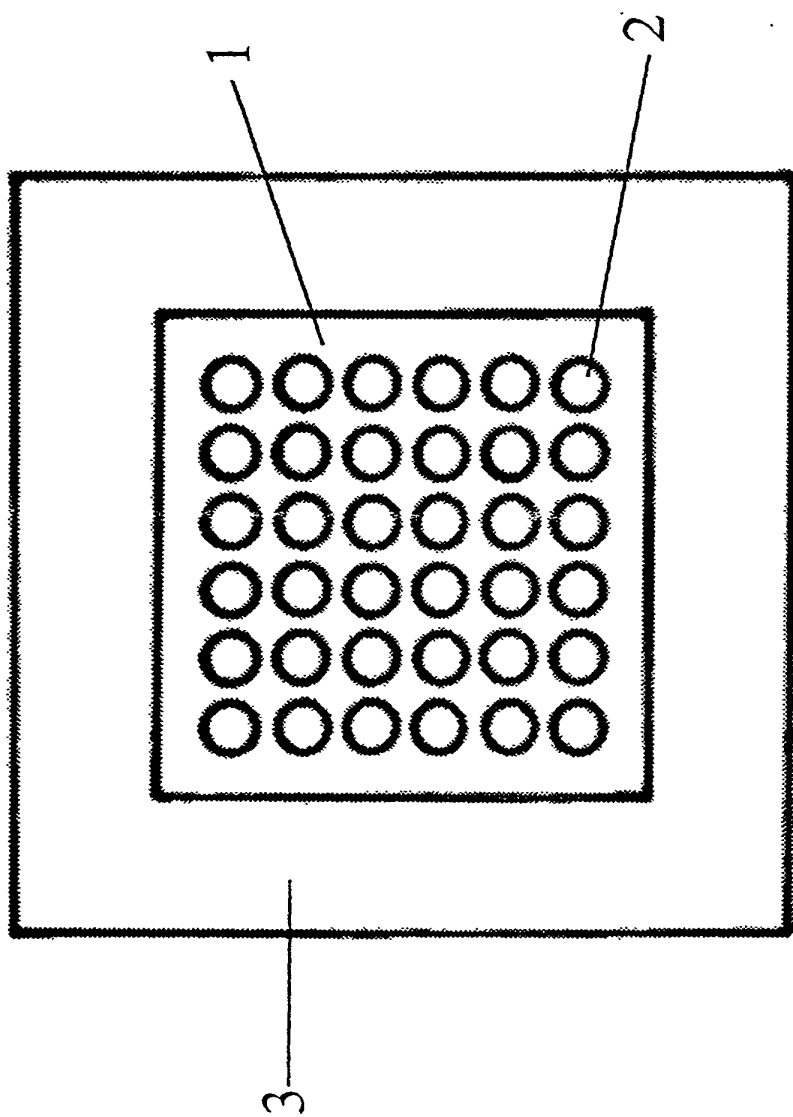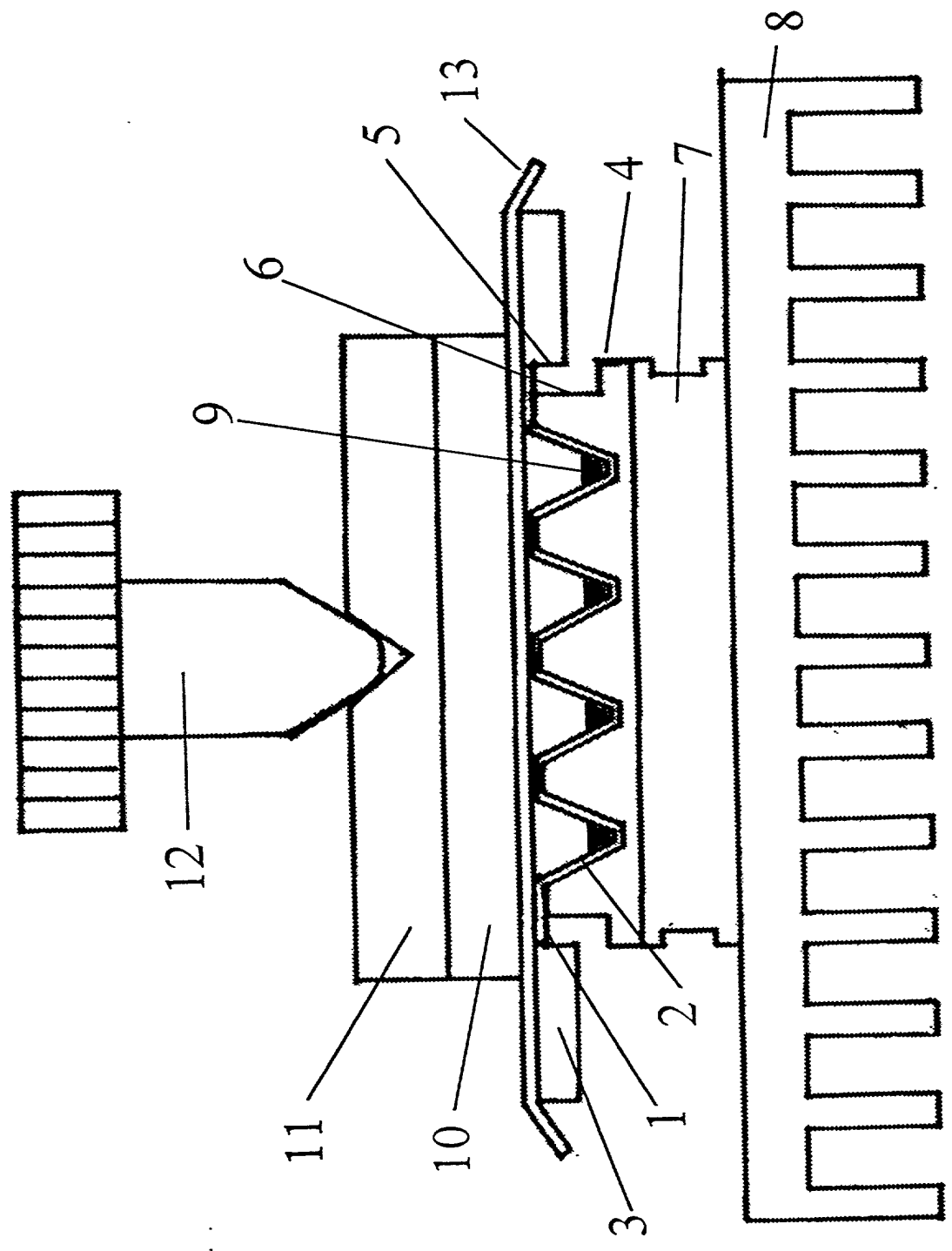

COMBINED DECLARATION FOR PATENT APPLICATION AND
POWER OF ATTORNEY

(Includes Reference to PCT International Applications)

Agency's Docket Number

As a below named inventor, I hereby declare that:

My residence, post office address and citizenship are as stated below next to my name,

I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter which is claimed and for which a patent is sought on the invention entitled:

ULTRATHIN-WALLED MULTIWELL PLATE FOR HEAT BLOCK THERMOCYCLING

the specification of which (check only one item below):

☐ is attached hereto.

☐ was filed as United States application

Serial No. _____

on _____,

and was amended

on _____ (if applicable).

☐ was filed as PCT international application

Number PCT/EP99/08178

on October 28, 1999

and was amended under PCT Article 19

on _____ (if applicable).

I hereby state that I have reviewed and understand the contents of the above-identified specification, including the claims, as amended by any amendment referred to above.

I acknowledge the duty to disclose information which is material to the examination of this application in accordance with Title 37, Code of Federal Regulations, §1.56(a).

I hereby claim foreign priority benefits under Title 35, United States Code, §119 of any foreign application(s) for patent or inventor's certificate or of any PCT international application(s) designating at least one country other than the United States of America listed below and have also identified below any foreign application(s) for patent or inventor's certificate or any PCT international application(s) designating at least one country other than the United States of America filed by me on the same subject matter having a filing date before that of the application(s) of which priority is claimed:

PRIOR FOREIGN/PCT APPLICATION(S) AND ANY PRIORITY CLAIMS UNDER 35 U.S.C. 119			
Country (if PCT indicate "PCT")	Application Number	Date of Filing (day, month, year)	Priority Claimed Under 35 USC 119
EP	98120187.4	29.10.1998	<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No
			<input type="checkbox"/> Yes <input type="checkbox"/> No
			<input type="checkbox"/> Yes <input type="checkbox"/> No
			<input type="checkbox"/> Yes <input type="checkbox"/> No
			<input type="checkbox"/> Yes <input type="checkbox"/> No

**COMBINED DECLARATION FOR PATENT APPLICATION AND
POWER OF ATTORNEY (Continued)**
(Includes Reference to PCT International Applications)

Attorney's Docket Number

I hereby claim the benefit under Title 35, United States Code, §120 of any United States application(s) or PCT international application(s) designating the United States of America that is/are listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in that/those prior application(s) in the manner provided by the first paragraph of Title 35, United States Code, §112, I acknowledge the duty to disclose material information as defined in Title 37, Code of Federal Regulations, §1.56(a) which occurred between the filing date of the prior application(s) and the national or PCT international filing date of this application:

PRIOR U.S. APPLICATIONS OR PCT INTERNATIONAL APPLICATIONS DESIGNATING THE U.S. FOR BENEFIT UNDER 35 U.S.C. 120:					
U.S. APPLICATIONS			STATUS (Check One)		
U.S. Application Number	U. S. Filing Date		Patented	Pending	Abandoned
PCT APPLICATIONS DESIGNATING THE U.S.					
PCT Application No.	PCT Filing Date	U.S. Serial Numbers Assigned (if any)			

POWER OF ATTORNEY: As a named inventor, I hereby appoint the following attorney(s) and/or agent(s) to prosecute this application and transact all business in the Patent and Trademark Office connected therewith.

Frank J. Jordan Reg. No. 20,456
C. Bruce Hamburg Reg. No. 22,389
Lainie E. Dolinger Reg. No. 36,123

Herbert F. Ruschmann Reg. No. 35,341
Marvin Turken Reg. No. 18,330
Alfred D'Andrea Reg. No. 27,752

Send Correspondence To:	Jordan and Hamburg 122 East 42nd Street New York, New York 10168	Direct Telephone Calls to: C. Bruce Hamburg (212) 986-2340
-------------------------	--	--

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under section 1001 of Title 18 of the United States Code, and that such willful false statements may jeopardize the validity of the application or any patent issuing thereon.

Full Name of Sole or First Inventor SALUZ, Hans-Peter	Inventor's Signature <i>[Signature]</i>	Date 12 April 2001
Residence Germany		Citizenship Swiss
Post Office Address Dorfstr. 22 D-07646 Oberbodnitz, DE DEX		

Full Name of Second Joint Inventor, if any TRETIKOV, Alexandre	Inventor's Signature <i>[Signature]</i>	Date 12 April 2001
Residence Germany		Citizenship Russia
Post Office Address Am Herrenberge 11, D-07745 Jena, DE DEX		